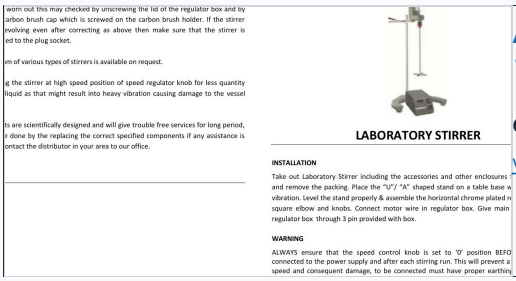
Importadora Gyb Limitada
Padre orellana 1476
Santiago CL13101 Chile

worm out this may checked by unscrewing the lid of the regulator box and by
arbon brush cap which is screwed on the carbon brush holder. If the stirrer
evolving even after correcting as above then make sure that the stirrer is
ed to the plug socket.

em of various types of stirrers is available on request.

g the stirrer at high speed position of speed regulator knob for less quantity
liquid as that might result into heavy vibration causing damage to the vessel

ts are scientifically designed and will give trouble free services for long period,
t done by the replacing the correct specified components if any assistance is
ontact the distributor in your area to our office.

LABORATORY STIRRER

INSTALLATION

Take out Laboratory Stirrer including the accessories and other enclosures
and remove the packing. Place the "U"/ "A" shaped stand on a table base w
vibration. Level the stand properly & assemble the horizontal chrome plated r
square elbow and knobs. Connect motor wire in regulator box. Give main
regulator box through 3 pin provided with box.

WARNING

ALWAYS ensure that the speed control knob is set to '0' position BEFO
connected to the power supply and after each stirring run. This will prevent a
speed and consequent damage, to be connected must have proper earthin

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