
STERLITECH

Corporation

Sepa CFX Cell - Titanium

Assembly & Operation Manual

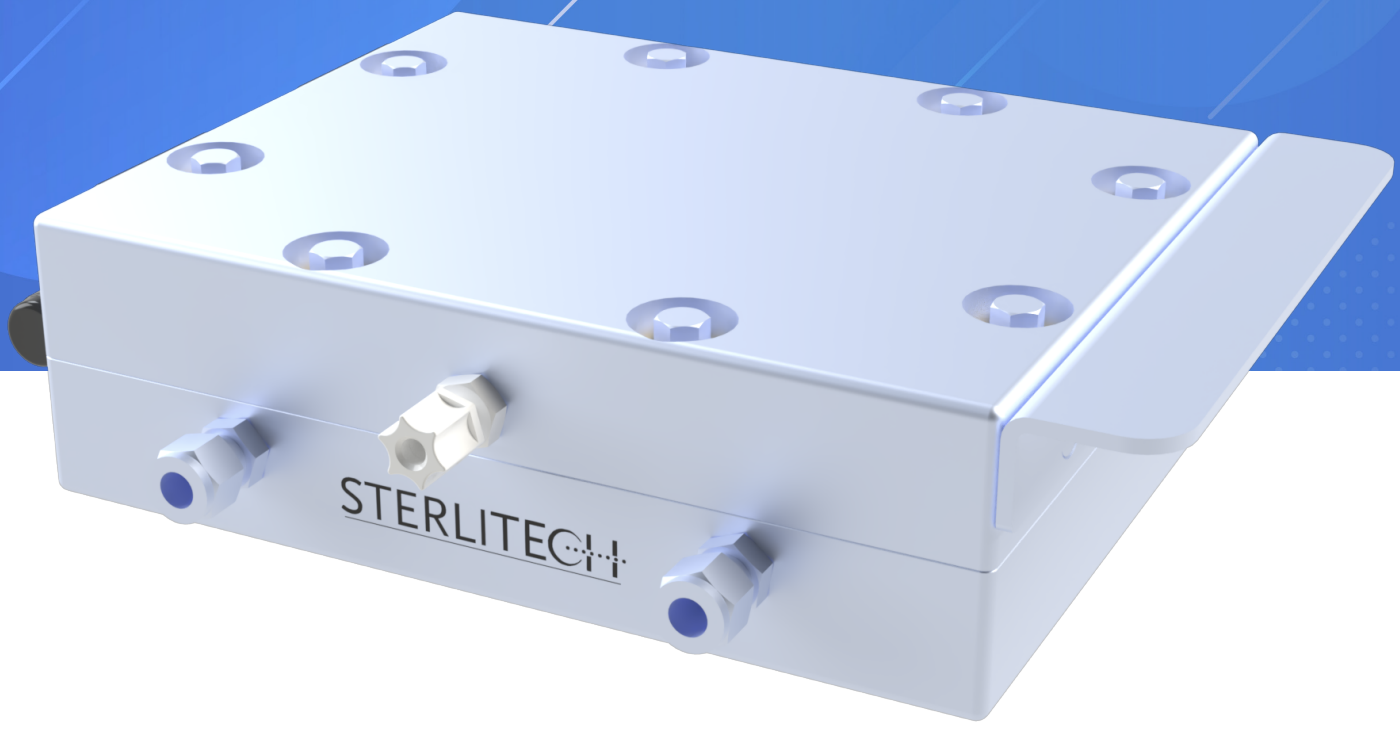

Copyright © 2026

Table of Contents

Introduction 3

Sepa CFX Cell Components 4

Sepa CFX Cell Assembly 5

Sepa CFX Cell System Assembly 6

Operation of the Sepa CFX Cell 7

Supplementary Operating Information 8

Accessory and Spare Part Ordering Information 9

Return Material Authorization 9

Warranty 9

Technical Assistance 9

 Copyright © 2026

Introduction

The Sepa CFX Cell is a laboratory-scale crossflow filtration unit designed to evaluate the performance of flat sheet membranes in a variety of applications. It simulates the flow dynamics of larger, commercially available membrane elements such as industrial spiral wound membrane elements. By using a combination shims, feed spacers, and membranes, users can vary the operating conditions and fluid dynamics over broad ranges.

Table 1 outlines the operational parameters and technical specifications of the Sepa CFX Cell.

Prior to operating or servicing this device, this manual must be read and understood.

Table 1: Sepa CFX Cell Features and Technical Specifications

PARAMETER	DESCRIPTION
Membrane Operation	Cross Flow
Cell Material	Grade 7 Titanium
Membrane Active Area	140 cm² (22 in²)
Maximum Pressure	139 bar (2000 PSI)
Maximum Temperature*	177°C (350°F)
Recommended Feed Rate	0.5-6 LPM (0.13-1.6 GPM)
Channel Dimensions	
Channel Depth	1.9 mm (0.075 in)
Channel Width	95.3 mm (3.75 in)
Fittings	
Inlet/Outlet Size (Tube OD)	3/8" Compression
Permeate Size (Tube OD)	1/4" Compression
Securing Method	High Pressure Screws


** Maximum Pressure and Temperature ratings are dependent on a variety of factors. Listed values are guidelines based upon material properties and cannot be guaranteed in every application.*

www.sterlitech.com info@sterlitech.com 1-877-544-4420
Copyright © 2026

Sepa CFX Cell Components

Figure 1 illustrates a typical assembly of the Sepa CFX Cell. The cell is also packaged with a length of permeate tubing and membrane cutting template.

Figure 1: Sepa CFX Cell Assembly

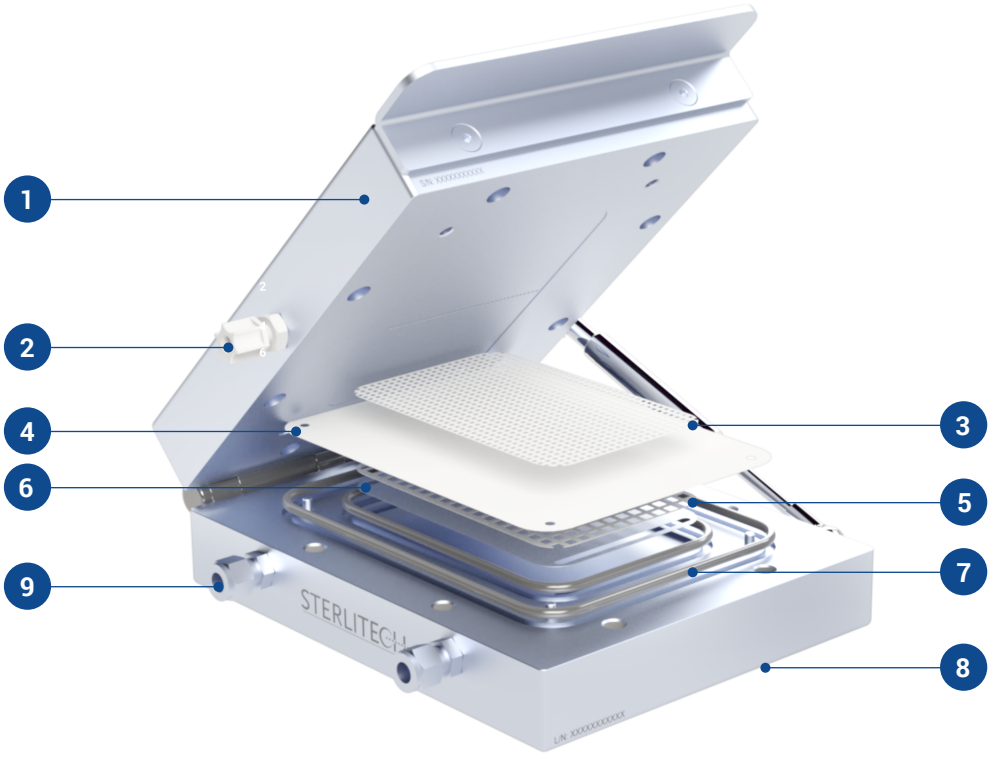

NUMBER	DESCRIPTION	MATERIAL	WETTED?
1	Cell Top	Titanium	Yes
2	Permeate Fitting	Polypropylene	Yes
3	Permeate Carrier	Polypropylene	Yes
4	Membrane	Various	Yes
5	Feed Spacer	Polypropylene	Yes
6	Shim **	316 Stainless Steel	Yes
7	O-Rings	PTFE Encapsulated Viton*	Yes
8	Cell Bottom	Titanium	Yes
9	Inlet/Outlet Fittings	Hastelloy C276	Yes
	Permeate Tubing	Polypropylene	Yes
	Screw & Washer Set	Stainless Steel / Brass	No

* *Standard, Other Materials Available* ** *Not Included, Available for purchase*

Copyright © 2026

Sepa CFX Cell Assembly

To assemble the Sepa Cell: (Video)

1. Wet the O-rings with a small amount of water or the fluid to be processed and install them into the grooves on the cell bottom. (Video)

   Caution: *Make sure that the O-rings lie flat in the grooves of the cell bottom. Leaking may occur if the O-rings do not lie flat in the grooves. The O-rings will be cut or crushed when the system is operating if the O-rings are not installed correctly.*

2. Optional: Install shim(s) into the cell bottom channel followed by the selected spacer. Spacers should be selected first, based on the desired thickness; shims should then be selected to fill the remaining channel depth. The overall channel depth is 0.075".

3. Place a piece of pre-cut membrane onto the cell bottom using the two alignment pins to hold the membrane in position. The membrane should be installed with the shiny side (active side) down, toward the cavity in the cell bottom.

   Note: *If you are cutting your own membranes, the membrane's outer edge should fit between the inner and outer O-rings in the cell body to prevent leakage.*

4. Wet the permeate carrier with water or the fluid to be processed and place it into the cavity in the cell top. The surface tension caused by wetting the carrier will keep it in place.

5. Gently close the cell top onto the cell bottom ensuring that the holes in the cell top align with the pins in the cell bottom.

6. Insert one of the brass washers into each of the eight holes on the cell top.

7. Secure the assembly by threading each of the provided screws from the top into the cell bottom. Tighten the screws in an alternating pattern to ensure even pressure on the O-rings.


   Note: *Pressure may need to be applied to the cell handle to keep the cell top from opening while securing with screws.*

www.sterlitech.com info@sterlitech.com 1-877-544-4420
Copyright © 2026

Sepa CFX Cell System Assembly

The Sepa Cell requires a cross flow system to operate, such as a Benchtop System or Skid System. Alternatively, **Figure 2** lists the key components required to operate a Sepa cell. These components are sold separately and can be found here.

Figure 2: Cross Flow System Diagram for Sepa CFX Cell

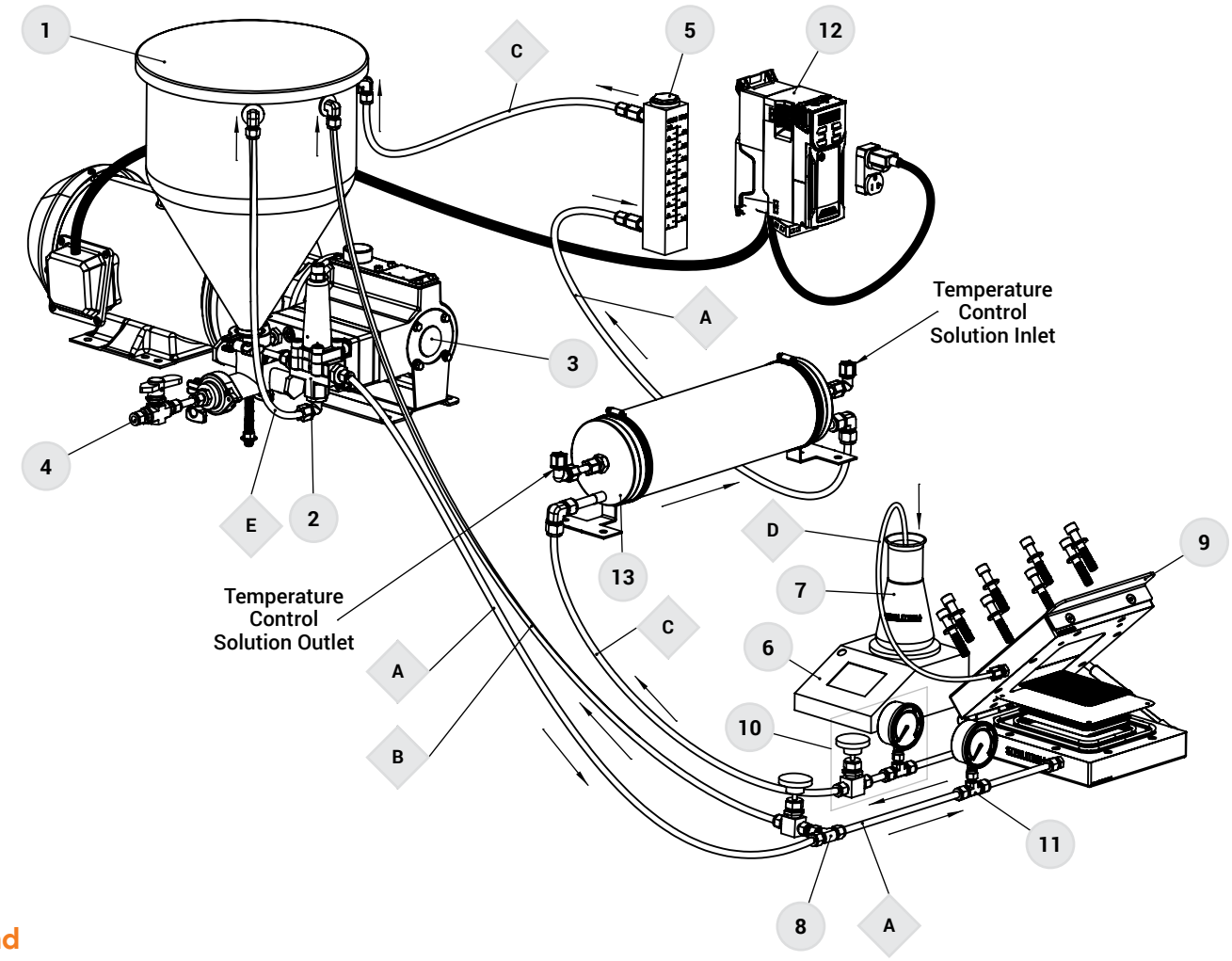

Legend

1	5 Gal Conical Feed Tank with Tee
2	Pressure Relief Valve
3	Motor & Pump Assembly
4	Drain Valve
5	Flow Meter
6	Scale For Measuring Permeate Flux
7	Permeate Collector
8	Bypass Needle Valve with Tee Tube Fitting
9	Sepa CFX Cell Assembly

10	Concentrate Control Valve Assembly
11	Feed Pressure Gauge
12	Variable Frequency Drive (VFD)
13	Heat Exchanger
A	High Pressure Flex Or Rigid Tubing, Feed Line
B	Low Pressure Tubing, Bypass Line
C	Low Pressure Tubing, Return Line
D	Low Pressure Tubing, Permeate Output Line
E	Low Pressure Tubing, Relief Return Line

Copyright © 2026

Operation of the Sepa CFX Cell

To operate the Sepa CFX Cell:

1. Open all the valves in the system to relieve any pressure build-up before turning on the feed pump.
2. Turn the system ON and set the feed pump to desired flow rate.

   Note: *If a bypass line is present, slowly close the bypass control valve; this valve is required only at low pressures and flows. Leaving the bypass line slightly open can also help to stabilize the system.*
3. Set the feed pressure by gradually closing the concentrate control valve assembly. Do not exceed the maximum pressure rating of the cell, holder, membrane or system.

   Note: *A good starting point may be to set the system at 1.7 bar (25 PSI) pressure for ultrafiltration.*
4. Verify the feed pressure reading on the concentrate pressure gauge.
5. Adjust the concentrate flow control valve to obtain the desired pressure and flow.

   Note: *Experimentation is required to determine the optimum settings of pressure, flow rate, and shim/spacer combination to use on the chosen membrane and the fluid being processed.*
6. Once the experiment is complete, shut down the system by turning the feed pump OFF and opening all the valves to relieve any pressure build up. Relieve the cell holder pressure by turning the relief valve knob counterclockwise.

To replace a membrane filter:


1. Turn the feed pump OFF.
2. Open all the valves in the system to relieve any pressure build up.
3. Unscrew all the screws from the cell top.
4. Open the cell by lifting the cell handle.
5. Remove the membrane.
6. Install the new membrane.
7. Close the cell top, checking that the cell top holes align with the cell bottom pins.
8. In an alternating pattern, secure the cell top to the cell bottom using all eight screws.

www.sterlitech.com info@sterlitech.com 1-877-544-4420
Copyright © 2026

Supplementary Operating Information

Additional Resources

- MPD Cell Quick Reference Guide
- MPD Video Channel
  - Cell Assembly
  - Shim and Spacer Installation
  - System Operation
- CF Filtration Handbook
  - Basics of Cross Flow Filtration
  - FAQs
  - Example Studies
  - Preconditioning Membranes
  - Temperature Effect on Membranes
- How to Swage a Compression Fitting
- How to Cut Custom Membrane Filters

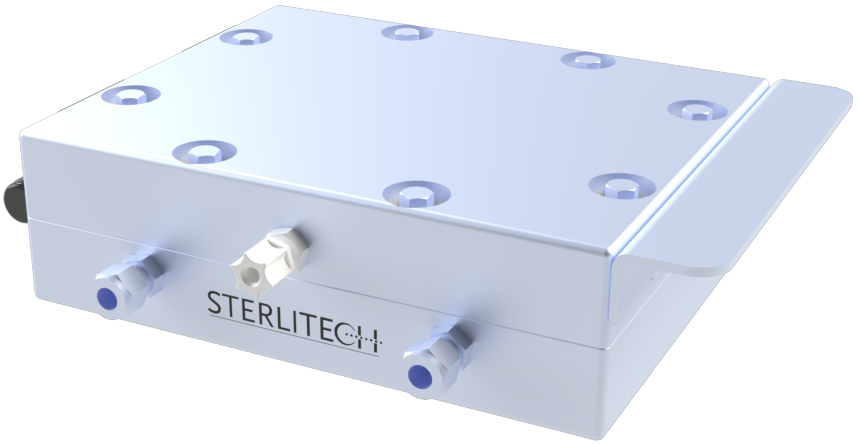

Recommended Blog Articles

- How to Calculate Cross Flow Velocity
- How to Select a Spacer
- How to Clean Membranes

Analyzing Membrane Performance Results

Spiral wound membranes containing mesh spacers are usually operated with a fluid velocity across the membrane surface of 0.1–0.5 m/sec (0.3–1.6 ft/sec). Higher velocities in membranes may lead to excessive pressure differential across the membrane and possible damage. Hollow fiber membranes, tubular membranes, and membranes with tubular spacers may be operated at higher velocities, but this may not contribute to a more effective operation.

Other parameters such as viscosity, pressure, and suspended solids may also affect performance/operation. Experimentation with the Sepa Cell can help predict the best operating parameters.

Cell and Membrane Cleaning

The Sepa Cell can be cleaned with or without the membrane installed. Leave the membrane installed to simulate the cleaning conditions of Clean-In-Place (CIP), a common practice in larger systems with spiral-wound or tubular membranes.

During CIP, cleaning solutions are re-circulated and, in some cases, allowed to sit for a period of time within the cell body. In some cases, the feed pump can be used to re- circulate the cleaning solutions.

Refer to the Membrane Cleaning Blog linked above for recommended cleaning instructions.

 Copyright © 2026

Accessory & Spare Part Ordering Information

PRODUCT	LINK
System Components (Pump, Flow Meters, Pressure Gauges, etc.)	System Components
Sepa Shim Assortment Pack	1231104
Sepa Feed Spacer Assortment Pack	1232558
Sepa Permeate Carrier Pack	1142817
Sepa Steel Rule Die	Sepa Die
Membrane Packs	Flat Sheet Membranes
Spare Parts (Alternative O-Rings, Sintered Support Plate, Knobs)	Spare Parts

Product Support

Return Material Authorization

If materials are to be returned to Sterlitech for repair, evaluation, or warranty consideration, a Return Material Authorization (RMA) number and form must be obtained from Sterlitech prior to the return. Contact Sterlitech's Customer Service Department for these forms.

Be sure to include a complete, detailed written reason for the return. Also, include serial numbers, installation and removal dates, and any other information that is available. The crossflow cells have a serial number imprinted on them.

Indicate the proposed disposition of the material, and reference the RMA number on all packages or cartons.All material must be" shipped to Sterlitech with freight prepared by the customer.

Warranty

The following is made in lieu of all other warranties expressed or implied. Sterlitech Corporation guarantees equipment to be free from defects in material and workmanship when operated in accordance with written instructions for a period of one year from receipt. Parts not manufactured by Sterlitech are covered by their manufacturer's warranties, which are normally for one year.

Manufacturers and Seller's only obligation shall be to issue credit against the purchase or replacement of equipment proved to be defective in material or workmanship. Neither manufacturer nor seller shall be liable for any injury, loss or damage, direct or indirect, special or consequential, arising out of the use of, misuse, or the inability to use such product.

The information contained herein is based on technical data and tests, which we believe to be reliable, and is intended for use by persons having technical skill at their discretion and risk. Since conditions of use are outside Sterlitech's control, we can assume no liability whatsoever for results obtained or damages incurred through the application of the data presented.

This information is not intended as a license to operate under, or a recommendation to infringe upon, any patent of Sterlitech or others covering any material or use.The foregoing may not be altered except by a written agreement signed by officers of the manufacturer.

Technical Assistance

Please contact us if you have any questions or technical inquiries about our products by calling Sterlitech Corporation at **1-877-544-4420** or by visiting www.sterlitech.com.

 Copyright © 2026

Company Background

Founded in 2001 in Kent, WA, Sterlitech Corporation manufactures and markets filtration-focused laboratory products to a broad spectrum of scientific and industrial sectors. Its line of flat sheet membranes and tangential flow cells deliver industry-leading performance reliable results. Configured for reverse osmosis, nanofiltration, ultrafiltration, and microfiltration applications, Sterlitech's benchtop cross flow membrane test system provides the versatility required to innovate.

Sterlitech's comprehensive line of products is supported by the expertise of its technical specialists who can assist with application-specific product selection, and provide customized solutions where necessary. Unique problem-solving approaches, flexibility, and consistent quality have made Sterlitech Corporation a renowned global provider of filtration products and equipment.

Sales Inquiries
sales@sterlitech.com

Accounts Payable
AP@sterlitech.com

Accounts Receivable
AR@sterlitech.com

Press Inquiries
PR@sterlitech.com



STERLITECH

4620 B Street NW, Suite 101
Auburn, WA 98001 USA
Tel: 877-544-4420 or 1-253-437-0844
Fax: 1-253-437-0845
Email: info@sterlitech.com

© 2026 Sterlitech Corporation